

Super FastPrep-2™



Owners Manual

Cat. No.
116012500

For Research Use Only

Protocol Revision #116012500-15OCT

mpbio.com/sampleprep

MP



™

Table Of Contents

Precautionary Instructions	3
Electrical Safety	3
Safety Equipment	3
Operation Precautions	3
Introduction	4
Product Contents	4
Picture of Contents	5
Super FastPrep-2™ Operating Instructions	6
Battery Charging	6
Holding the <i>Super FastPrep-2™</i>	6
Sample Processing	7
Cleaning	9
Warranty and Disclaimer	10

It is recommended to keep the original packaging in case you need to return the product for any reason to your distributor or MP Biomedicals.

Note: MP Biomedicals is not responsible for any harm or damage resulting directly/indirectly from the use of the *Super FastPrep-2* or related product components, materials, or technologies.

Some reagents in kits are irritants. Wear protective gloves and eye protection and follow the safety guidelines and rules enacted by your research institution or facility. Refer to product/reagent specific MSDS information from the manufacturer.

PRECAUTIONARY INSTRUCTIONS

Be sure to read all safety instructions in this manual carefully and thoroughly. Pay attention to warnings before using the *Super FastPrep-2*. Failure to follow all instructions may lead to injury.

⚠ Warning: Indicates a potential hazard that should be avoided.

Electrical Safety

- ⚠ Avoid contact with grounded surfaces.
 - ⚠ Outlets for charging must match the plug on the battery charging station.
 - ⚠ A voltage converter is included for countries that use 220-240 VAC.
- ⚠ Do not use the *Super FastPrep-2* in the rain or under wet conditions.
- ⚠ Do not misuse the cord.

Safety Equipment

- ✓ The *Super FastPrep-2* can be very loud so it is recommended to use ear protection during operation.
- ✓ Always use eye protection while operating the *Super FastPrep-2*.
- ⚠ For further safety information regarding the rotary tool, battery, and charging station please see accompanying literature.

Operation Precautions

- ⚠ Operation of the *Super FastPrep-2* with a plastic tube at full speed for more than 10 seconds may cause the grinding matrices within the tube to overheat, leading to sample degradation and/or tube failure. If grinding of a sample for over 10 seconds is desired, lower speeds should be used and/or periodic stopping and cooling of the tube should be performed.
- ⚠ When using metal or ceramic grinding matrices greater than 2mm in diameter, do not operate at a speed greater than 25 (4,000 cycles per minute) or tube failure may result.
- ⚠ Although the modular design of the *Super FastPrep-2* allows it fit on many off-the-shelf rotary tools from different manufacturers, using the *Super FastPrep-2* with a rotary tool with an electric stop (stops immediately after switching it off) will damage the *Super FastPrep-2* and void the warranty. Use the *Super FastPrep-2* only with the rotary tool supplied, which has been specifically modified to disable its electric stop and coast to a stop when switched off.

INTRODUCTION

The *Super FastPrep-2* is an attachment for a rotary tool developed for vigorous cell disruption (bead beating) and allows the researcher/investigator to process samples in the field, at the site of sample collection. It can be used to lyse microbes in soil, sediment, sludge, and fecal samples and can effectively process tough-to-lyse fungal, algal, plant, and animal tissues. It can be used at any remote location and in most weather conditions when immediate sample collection, processing, and preservation are required by the researcher. The device is compatible with most 2 mL tubes containing lysis matrix.

The *Super FastPrep-2* consists of a unique, reciprocating grinding mechanism mounted on a rotary tool. It is designed to operate at 600 to 4,800 cycles per minute as long as charge remains in the battery. This allows the researcher to maintain control of the speed and time of sample lysis. Overall, the *Super FastPrep-2* is ideal for cordless, portable, and convenient processing of tough-to-lyse biological and environmental samples, both in the lab and in the field.

PRODUCT CONTENTS

- 1) One (1) *Super FastPrep-2*
- 2) One (1) Rotary tool
- 3) Two (2) Lithium Ion Rechargeable Batteries
- 4) One (1) Battery Charging Station
- 5) Ten (10) sample tubes
- 6) Four (4) tube support plugs
- 7) Product Manuals
- 8) 220-240 VAC to 100-120 VAC Voltage Converter to countries where necessary

Super FastPrep-2™



OPERATING INSTRUCTIONS

Battery Charging

Connect the battery charging station to an appropriate power outlet (use the power converter to use the charging station outside of North America)

Place the battery in the charging station by aligning the raised ribs of the battery pack with the slots in the charging station. The light on the charging station should flash green while the battery is charging. The battery pack *is fully charged when the green light stops flashing and remains on.*

Holding the *Super FastPrep-2*

During operation, hold the *Super FastPrep-2* as in the picture below.



The angular position of the rotary tool relative to the *Super FastPrep-2* may be adjusted by removing the 3 screws at the base of the *Super FastPrep-2*, rotating the adapter to the desired position, and then replacing the 3 screws. See the picture below.



Sample Processing


✓ **Sample Sources** – The *Super FastPrep-2* device can be used to effectively lyse microbes in soil, sediment, sludge, and fecal samples and can effectively process tough-to-lyse fungal, algal, plant, and animal tissue. The *Super FastPrep-2* is compatible with skirted, 2 mL tubes. Two mL tubes having conical bases are compatible with the *Super FastPrep-1*, however support plug (included) must be installed to help prevent failure of those tubes.

✓ **Operating Temperatures** – (Optimum -20°C to 45°C) Severe cold conditions can drain the battery and affect the number of samples that can be processed with the *Super FastPrep-2*. Keeping the battery warm (e.g., in the pocket of a jacket or parka) prior to use can extend battery life in colder climates and operating conditions.

✓ **Processing Time and Speed** – (Generally <5 seconds at 4,000 cycles per minute) The extent of sample lysis will vary according to the processing speed and time. Due to the extreme pulverizing forces generated by the *Super FastPrep-2*, the temperature of samples can increase as a function of lysis speed and time. This may

affect the integrity of DNA/RNA/protein released from biological materials. If longer processing times (>5 seconds) are required, it is recommended to let the sample cool on ice for a minimum of two minutes before continuing with sample processing and/or to reduce the lysis speed.

1. Remove *Super FastPrep-2* and rotary tool from the box or carrying case.
2. Screw rotary tool into *Super FastPrep-2* to a tight fit.
3. Open lid and insert a screw-cap 2mL tube into the holder making sure that it snaps in securely. If you are using one or more 1/4" diameter metal or ceramic grinding balls with a conical-bottomed plastic tube, use the support plug to help prevent failure of the tube. Hold the tube and support plug between thumb and index finger. Press them together into the holder of the *Super FastPrep-2*.

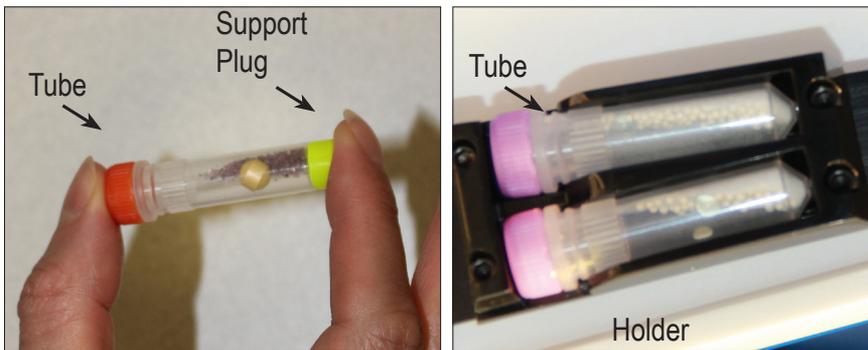

- ⚠ It is recommended that for best results, only Lysing Matrix Tubes supplied by MP Biomedicals be used in the *Super FastPrep-2*.
- ⚠ If you use a 2 mm (or greater) diameter metal or ceramic ball in a plastic tube, do not run the rotary tool above the speed setting of 25 (about 4,000 cycles per minute). Doing so may cause failure.
- ⚠ If you use a tube with a tethered cap, be sure to insert the tube into the holder with the plastic loop facing down to prevent it from damaging the clear window on the lid to the *Super FastPrep-2* during operation.

4. Close the lid and set the speed control on the rotary tool to the desired speed. Below is a table of the approximate *Super FastPrep-2* speed (+/- 1,000 cycles per minute) for each corresponding position of the speed control slide on the Dremel® 8220 rotary tool, running with an empty tube.

Slide	Cycles Per Minute
5	600
10	1,400
15	1,900
20	2,300
25	4,000
30	4,800

5. While holding the lid closed with one hand, flip the on/off switch on the rotary tool with the other hand to run the *Super FastPrep-2* for the desired period of time.

⚠ Caution: Do not operate the *Super FastPrep-2* while the lid is open.

6. After turning off the rotary tool, open the lid and remove the 2mL tube. Continue with downstream purification/isolation procedure.

Refer to corresponding purification kit instructions for details concerning subsequent isolation of DNA/RNA/Protein from lysate material. For a complete listing of FastDNA, FastRNA, and FastProtein purification and isolation kits please visit www.mpbio.com/sampleprep.

Cleaning the *Super FastPrep-2*

⚠ Make sure to disconnect the rotary tool from the *Super FastPrep-2* before cleaning. Wear protective gloves (latex/nitrile) during all cleaning procedures.

Clean the *Super FastPrep-2* as necessary. The exterior of the *Super FastPrep-2* can be cleaned with a damp (water, dilute ethanol) paper towel or similar. Clean the lysis chamber (exposed when lid is open) by wiping down the inner walls and tube holder. Allow all components to air dry prior to use.

⚠ Should a lysis tube fracture or break in the sample chamber during processing, stop processing and immediately clean it as recommended above.

For Research Use Only

All of our products are intended for research use only. Not for use in diagnostic procedures.

Warranty

This product is guaranteed for one (1) year from the time of receipt. See below for repair/return policy.

Returns

For repairs made under warranty, the Goods must be returned to MP Biomedicals intact (non-tampered with, non-disassembled) with all components together and in their original packaging. Goods may not be returned for credit except with Seller's (i.e., MP Biomedicals) permission, and then only in strict compliance with Seller's return shipment instructions. Returned Goods (other than defective products) must be returned freight prepaid by the Customer and will not be accepted without prior authorization of MP Biomedicals. The Goods must be returned in their original packaging and in resalable condition. Any returned items may be subject to a 20% processing (restocking) fee.

Terms of Purchase

We warrant to you, our direct customer, that our goods shall conform substantially to the description of such goods as provided in our catalogs and literature accompanying the goods until their respective expiration dates. THIS WARRANTY IS EXCLUSIVE, AND WE MAKE NO OTHER WARRANTY, EXPRESS OR IMPLIED, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE. Our warranty shall not be effective if we determine, in our sole discretion, that you have altered or misused the goods or have failed to use or store them in accordance with instructions furnished by us. Our sole and exclusive liability and your exclusive remedy with respect to goods proved to our satisfaction (applying analytical methods reasonably selected by us) to be defective or nonconforming shall be the replacement of such goods free of charge, upon the return of such goods in accordance with our instructions, although at our discretion we may provide a credit or refund. IN NO EVENT SHALL WE BE LIABLE UNDER ANY LEGAL THEORY (INCLUDING BUT NOT LIMITED TO CONTRACT, NEGLIGENCE, STRICT LIABILITY IN TORT OR WARRANTY OF ANY KIND) FOR ANY INDIRECT, SPECIAL, INCIDENTAL, CONSEQUENTIAL OR EXEMPLARY DAMAGES (INCLUDING BUT NOT LIMITED TO LOST PROFITS), EVEN IF WE HAD NOTICE OF THE POSSIBILITY OF SUCH DAMAGES. If we manufacture custom goods for you based on instructions, specifications, or other directions you provide to us, we shall not be liable for the lack of sufficiency, fitness for purpose or quality of the goods to the extent attributable to such instructions, specifications, or other directions. We shall not be liable for any loss, damage or penalty as a result of any delay in or failure to manufacture, deliver or otherwise perform hereunder due to any cause beyond our reasonable control.

Patent/Trademark Disclaimer

Some technologies of the *Super FastPrep-2* sample processor are patent pending. Unless stated otherwise, no license or immunity under any patents shall be either granted or implied by the sale of MP Biomedicals products. Dremel® is a trademark of Robert Bosch Tool Corporation.

Protocol Revision # 116012500-15OCT

Cat. No. 116012500

Worldwide Ordering and Technical Support

USA

Worldwide Headquarters

Tel: 800.854.0530
Fax: 800.334.6999

Europe

Phone: 00800.7777.9999
Fax: 00800.6666.8888

Australia

Tel: +61 2.9838.7422
Tel: 1800.249.998
Fax: +61 2.9838.7390
Fax: 1800.818.671

Belgium

Tel: +32 2.466.00.00
Fax: +32 2.466.26.42

Canada

Tel: 800.854.0530
Fax: 800.334.6999

China

Tel: +86 21.6281.5005
Fax: +86 21.6281.5663

France

Tel: +33 3.88.67.54.25
Fax: +33 3.88.67.19.45

Germany

Tel: 0800.426.67.337
Fax: 0800.629.67.337

India

Tel: +91 22.22631257/58
Fax: +91 22.22631259

Japan

Tel: +81 03.3808.2102
Tel: 0120.78.8020
Fax: +81 03.3808.2401

Latin America

Tel: +1 440.337.1200
Fax: +1 440.337.1180

The Netherlands

Tel: 0800.0227.416
Fax: 0800.0227.489

New Zealand

Tel: +64 9.912.2460
Fax: +64 9.838.4209

Poland

Tel: +48 22.659.58.95
Fax: +48 22.658.45.05

Russia

Tel: +7 495.661.0008
Fax: +33 3.88.67.19.45

Serbia

Tel: +381 11.242.1972
Fax: +381 11.242.1949

Singapore

Tel: +65 6775.0008
Fax: +65 6774.6146

Switzerland

Tel: +41 61.271.0007
Fax: +41 61.271.0084

United Kingdom

Tel: 0800.282.474
Fax: 0800.614.735

mpbio.com/sampleprep
support@mpbio.com

MP Biomedicals • 29525 Fountain Parkway • Solon, OH 44139 USA • tel: 800.854.0530 • fax: 800.334.6999

MP



™